

interscience



MiniMix® 100 CC®

Lab Blender

JumboMix® Brochure N° 30910-0

GB

Simple, Accurate, Space-saving

The most compact blender in the world, BagMixer® 100 MiniMix® is now available in Click & Clean® version!

MiniMix® 100 automatically blends any solid sample in a sterile bag and allows safe and easy extraction of bacteria. You can pipet and analyse the filtrate without any risk of alteration or cross contamination. **This is the safest and most economical sample preparation method.**

The new MiniMix® 100 CC® incorporates the Click & Clean® and side-by-side paddle stopping systems, already highly acclaimed on BagMixer® 400 CC® and the adjustable blending strength system.



1 Multifunction screen

For easy control of the blending parameters

Blending time


Blending speed



2 Q-Tight® system


A special closing system that avoids every risk of leaking.

Leak-tight

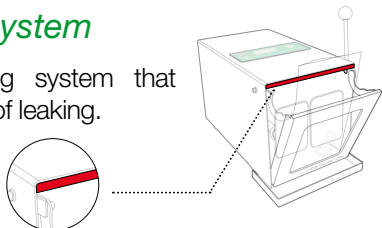

3 Adjustable blending strength PATENTED

You can adapt the blending chamber to the size of the sample to get OPTIMAL EXTRACTION!

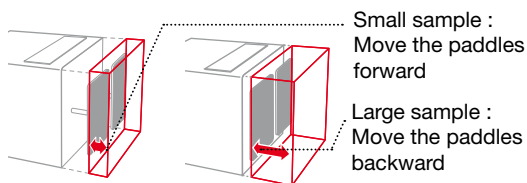


Small sample :
Move the paddles forward


Large sample :
Move the paddles backward

4 Click & Clean® system PATENTED


UNIQUE ! TOTAL CLEANING of the blending chamber is now possible !

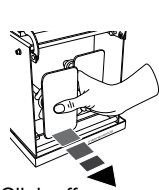

Click off

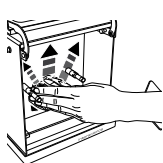


Clean

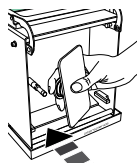

Click back on

5 Side-by-side paddle stop PATENTED

For easy and fast insertion of your sterile bags in the blender.

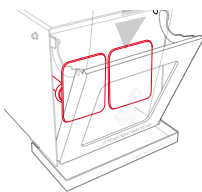

Advantages

- Compact, silent
- User friendly
- Optimal extraction
- Lifetime warranty on window door and shock absorbers
- Stainless steel
- Q-Tight® : Leak-tight
- Easy and total cleaning

Applications

- Food microbiology
- Pharmaceutical
- Veterinary, clinical (biopsies, organ sampling...)
- Molecular extraction (proteins, lipids, DNA...)
- Toxin research (HPLC, GCMS...)

Specifications

- Useful capacity : 5 - 100 ml
- For samples : 0.5 g - 10 g
- Adjustable time : 1 sec. à 59 min. countdown & open running (∞)
- Speed calibrated by microprocessor (3 - 12 strokes/sec.)
- Side by side paddle stop
- Q-Tight® seal : thickness of 12mm
- Stainless steel
- Drip tray
- Adjustable blending strength
- Click & Clean® system : removable & autoclavable paddles
- Multifunction screen
- Lifetime warranty on window door & shock absorbers

- Pwr. : 220 / 100 - 110 V – 50 / 60 Hz
- Dim. (LxWxH) : 37 x 20 x 20 cm
- Net weight : 13 Kg

Ref : 013 230



WEEE Compliant
Directive 2002/96/EC

RoHS Compliant
Directive 2002/95/EC

Please return this form to your local distributor or to interscience.



Also available

MiniMix® 100 ml blender

- Quotation
- Demonstration

I would like to receive more info on:

- BagMixer® 400 ml blender
- JumboMix® 3500 ml blender
- BagFilter® filter bags
- BagTools® accessories
- Gravimat® gravimetric dilutor
- easySpiral® automatic plater
- Scan® colony counters

Name.....

Company.....

Address.....

Zip/City..... Country.....

Phone..... Fax.....

email.....

web.....

interscience

Worldwide: 30, Ch Bois Arpents, 78860 St Nom FRANCE
Tel. +33(0)1 34 62 62 61 / Fax. +33(0)1 34 62 43 03
info@interscience.fr / www.interscience.fr

USA/Canada: 199 Weymouth St. ROCKLAND, MA 02370 USA
Tel. +1 781 792 2133 / Fax. +1 781 792 2134
info@intersciencelab.com / www.intersciencelab.com



MiniMix® 100 P CC®

- Variable blending speed (3-12 strokes/sec)
- Variable blending time (1 - 1h or ∞)
- Available in CoolMix® version


Get the best results with



Gravima®

Gravimetric dilutor

Automatically dilutes with the appropriate mass of solvent.

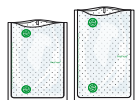

BagPage® 80 & 100

Full page sterile filter bags


Robust & food compatible, multilayer.



BagLight® 100

Plain sterile bags

Made of PolySilk®, extremely robust. For samples that do not require any filtration.

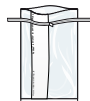

RollBag® 450

Sampling bag

For collecting, sampling & preservation of any sample. Stands up alone.



BagOpen® 100

Bag opener

Easy opening, extremely stable. Resistant to cleaning agents.



BagRack® 100

Storage rack

Stainless steel, autoclavable.



BagClip® 100

Closing clips

100% watertight & waterproof, closes any kind of bag.

MiniMix® Brochure N° 30010-0 - MiniMix® 100 CC®, Click & Clean® system, MiniMix® 100 P CC®, Gravimat®, BagPage®, BagLight®, RollBag®, BagOpen®, BagRack®, BagClip®, BagMixer®, JumboMix®, BagFilter®, BagTools®, easySpiral® and Scan® are registered trademarks of Interscience. Information in this leaflet is not contractually binding.